

## **VERSATILITY & VISIBILITY**

The Reform® Modular Pedicle Screw System provides surgeons increased flexibility, versatility and visibility to meet the varying requirements of degenerative and trauma procedures. Reform features a cobalt chrome tulip, a titanium triple lead proximally tapered thread to deliver strength, stability and efficiency to all thoracolumbar constructs. Modular tulips are available in standard and reduction configurations. Testing shows that the Reform Modular intraoperatively assembled construct has a pull off strength as strong as a Reform preassembled screw providing optimal security.







## REFORM® MODULAR PEDICLE SCREW

**Modular Screw Driver** offers a secure. low profile interface optimizing maneuverability, control, and dependability.


**Semi Aggressive Bone Planar** approximates the outer diameter of the tulip to allow for bone decortication to promote fusion and proper seating of the tulip

Low Profile Tulip Clamp provides a stable and secure interface and is available in straight and angled configurations. An audible click ensures that the tulip is securely assembled

**Polyaxial Screws** provide 30° angulation in all planes for intraoperative flexibility.

**Proximal Tapered Thread** design increases Bone Screw interface, enhancing pull-out strength.

**Triple Lead Thread** allows for efficient screw delivery.

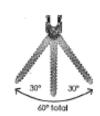









Bone Planar



Tulip Clamp

**Cobalt Chrome Tulip** reduces risk of head splay and increases overall construct stabilization. Standard and reduction tulip options are available.

**Square Threaded Locking Cap** geometry reduces risk of cross threading









## Reform Modular Sizing

| 4.5mm    | 5.5mm    | 6.5mm    | 7.5mm    | 8.5mm    | 9.5mm    |
|----------|----------|----------|----------|----------|----------|
| Diameter | Diameter | Diameter | Diameter | Diameter | Diameter |
| 25mm     | 30mm     | 30mm     | 35mm     | 40mm     | 60mm     |
| 30mm     | 35mm     | 35mm     | 40mm     | 45mm     | 70mm     |
| 35mm     | 40mm     | 40mm     | 45mm     | 50mm     | 80mm     |
| 40mm     | 45mm     | 45mm     | 50mm     | 55mm     |          |
| 45mm     | 50mm     | 50mm     | 55mm     | 60mm     |          |
|          |          | 55mm     |          | 70mm     |          |
|          |          |          |          | 80mm     |          |

## Reform Rod Sizing - 5.5mm ø

| Ti       | Ti              | CoCr     | CoCr            |
|----------|-----------------|----------|-----------------|
| Straight | Curved          | Straight | Curved          |
| 80mm     | 35-80mm (5mm)   | 80mm     | 35-80mm (5mm)   |
| 100mm    | 90-120mm (10mm) | 100mm    | 90-120mm (10mm) |
| 120mm    |                 | 120mm    |                 |
| 200mm    |                 | 200mm    |                 |
| 300mm    |                 | 300mm    |                 |
| 400mm    |                 | 400mm    |                 |
| 500mm    |                 | 500mm    |                 |

Precision Spine, Inc. 2050 Executive Drive, Pearl, MS 39208 Customer Service: 1.888.241.4773 Phone: 601.420.4244 Toll Free: 877.780.4370 Fax: 601.420.5501 www.precisionspineinc.com

